


waterpik® waterflosser®

QUICK START GUIDE WF-02


Visit
Waterpik.com/WelcomeWF
to watch How To Video
Register your product and
receive over \$50 in savings.
Waterpik.com/registration

1 Unit Prep (First time use only)


1. Remove reservoir

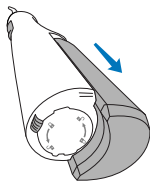


2. To remove battery door, twist, and pull.

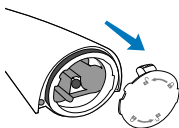


3. Insert 3 AA batteries.

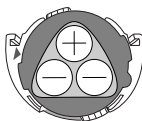


4. To replace battery door, line-up arrows, push and twist.

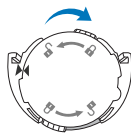

5. Insert tip firmly until it clicks into place.

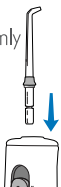

6. Fill reservoir with warm water.

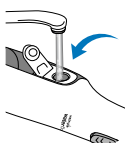

Important: Use Petroleum jelly to lubricate the battery cap sealing surface at every battery change. Apply a thin even film of Petroleum Jelly to the battery cap o-ring and the battery compartment sealing surface.

2 Water Floss




Slide power switch to low (I) position. Start with lowest setting and increase to your liking.



Close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.

Technique

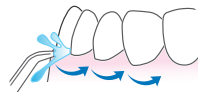

Aim water at the gumline at a 90 degree angle.

Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.


waterpik® waterflosser®

GUÍA RÁPIDA PARA COMENZAR WF-02


Visite
Waterpik.com/Welcome
para ver el video práctico
Registre su producto y reciba
más de USD 50 en descuentos.
Waterpik.com/registration

1 Preparación de la unidad

1. Quite el depósito de agua.

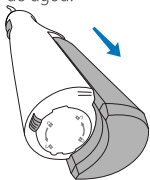


2. Para quitar la puerta de la batería, gírela y jale.

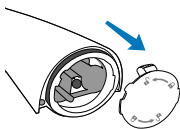

3. Introduzca 3 pilas AA




4. Para volver a colocar la tapa de las pilas, alinee las flechas, empuje y gire.

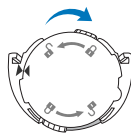

5. Inserte el cabezal con firmeza hasta que encaje en su lugar.




6. Llene el depósito con agua tibia.

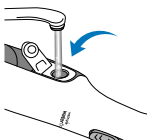

2 Cómo utilizar el Irrigador oral




Deslice el interruptor de encendido hacia la posición baja (I). Inicie con el ajuste más bajo e increméntelo según su preferencia.

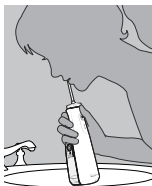

Cierre los labios lo suficiente para evitar salpicar mientras el agua continúa fluyendo de la boca hacia el lavabo.

Técnica

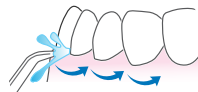

Dirija el agua hacia las encías en un ángulo de 90 grados. Siga las encías y haga pausas breves entre los dientes.

**EVITE DERRAMES: COLOQUE EL CABEZAL EN LA BOCA
ANTES DE ENCENDER LA UNIDAD**